

NUCLEIC ACID PURIFICATION
PURE AND SIMPLE®



Ionic® Purification System

SITE PREP GUIDE

REV 4.0

PART #: PUR-DOC-1

Purigen Biosystems, Inc.

5700 Stoneridge Drive, Suite 100
Pleasanton, CA 94588

TEL: +1 877 787-4436 **EMAIL:** info@purigenbio.com



Table of Content

Copyright.....

Support

Introduction

Additional Resources.....

Delivery and Installation

Box Dimensions and Contents.....

Laboratory Requirements.....

General Requirements.....

Instrument Dimensions.....

Placement Requirements.....

Lab Bench Guidelines.....

Vibration Guidelines

Electrical Requirements.....

Power Specifications.....

Recommendations

Power Cords.....

Fuses

Uninterruptible Power Supply

Product Certifications and Compliance

European CE Standards

Environmental Considerations

Pre-installation Checklist

Acknowledgment and Authorization

3

3

4

4

4

5

5

5

6

7

7

8

8

8

9

9

9

9

9

10

10

11

12

Copyright

© 2021 Purigen Biosystems, Inc. All rights reserved. Duplication and/or reproduction of all or any portion of this document without the express written consent of Purigen Biosystems, Inc. is strictly forbidden. Nothing contained herein shall constitute any warranty, express or implied, as to the performance of any products described herein. Any and all warranties applicable to any products are set forth in the applicable terms and conditions of sale accompanying the purchase of such product. Purigen Biosystems provides no warranty and hereby disclaims any and all warranties as to the use of any third-party products or protocols described herein. The use of products described herein is subject to certain restrictions as set forth in the applicable terms and conditions of sale accompanying the purchase of such product. “Ionic®” is a registered trademark of Purigen Biosystems, Inc. All other trademarks are the property of their respective owners. All products and services described herein are intended **FOR RESEARCH USE ONLY and NOT FOR USE IN DIAGNOSTIC PROCEDURES.**

All product names, brands, and trademarks are properties of their respective owners.

Support

For technical information or advice, please contact Purigen Biosystems Support at any time.

Phone: 1-877-PURIGEN (787-4436)

Email: support@purigenbio.com

Purigen Biosystems, Inc.

5700 Stoneridge Drive, Suite 100

Pleasanton, CA 94588

USA

Introduction

This guide contains specifications and guidelines for preparing your site for installation and operation of the Ionic® Purification System. Prior to installation of your Ionic system, a Purigen Biosystems Field Service Engineer (FSE) will contact you to schedule the install and confirm that your site meets the preparation requirements described in this document. If your system installation has already been scheduled, and you believe your site may not be ready in time, please contact your Purigen Biosystems Field Service Engineer (FSE) for assistance.

- Laboratory space requirements
- Electrical requirements
- Environmental constraints
- Computing requirements
- User-supplied consumables and equipment

Additional Resources

The following documentation is available for download from the Purigen Biosystems website.

TABLE 1: Additional documentation resources

Guide Name	Description
Ionic® Purification System User Guide	Provides an overview of instrument components and software, and procedures for proper instrument maintenance and troubleshooting
Ionic® Cells to Pure DNA Kit	Provides steps to use the Protocol and instructions for the Ionic Cells to Pure DNA Kit
Ionic® FFPE to Pure DNA Kit	Provides steps to use the Protocol and instructions for Ionic FFPE to Pure DNA Kit

Delivery and Installation

Following your purchase of an Ionic system, a Purigen Biosystems representative will contact you to confirm receipt of your order and provide an early estimate of delivery. Prior to delivery, a Purigen Biosystems FSE will contact you to schedule the installation.

The Ionic system will be shipped to you directly. In this case, you will receive shipment tracking information from Purigen Biosystems. Inspect the packages containing your Ionic system and report any signs of damage or rough handling to your Purigen Biosystems Support.



IMPORTANT

Do not unpack the instrument prior to installation. The system must be unpacked by a Purigen Biosystems FSE during installation.



CAUTION

Only Purigen Biosystems personnel can unpack, install, or move the instrument. Mishandling of the instrument can affect the alignment or damage instrument components.



CAUTION
After your Purigen Biosystems representative has installed and prepared the instrument, do not relocate the instrument. Moving the instrument improperly can affect the alignment and compromise yield and purity. If you want to relocate the instrument, contact your Purigen representative.

Box Dimensions and Contents

The Ionic system is shipped in a double-wall corrugated fiberboard box with a designer tray with full overlap, poly straps and hand-holds. Use the following dimensions to determine the minimum door width required to accommodate the shipping container.

TABLE 2: Box dimensions

Specification	Value
Height	62.7 cm (24.7 in.)
Width	47.2 cm (18.6 in.)
Depth	55.1 cm (21.7 in.)
Weight	29.0 kg (64 lbs.)

The box contains the instrument along with the following components:

- 1 power cable

Laboratory Requirements

This section provides requirements and guidelines to set up your lab space for the Ionic system properly. For more information, see ["Environmental Considerations"](#) on page 10.

General Requirements

- Adequate floor space near the installation site to accommodate the unpacking and unboxing of system components.
- Adequate clearance around the instrument for operator and service access.
- Access to receptacle for power.
- Rigid table or benchtop with the appropriate footprint and construction required to minimize unwanted vibration or instability.
- Adequate ventilation to accommodate heat output for all devices in the laboratory where system will be used.
- Adequate management of the laboratory environment to minimize fluctuations in temperature and humidity.

Instrument Dimensions

FIGURE 1: Instrument dimensions with leveling system (front)


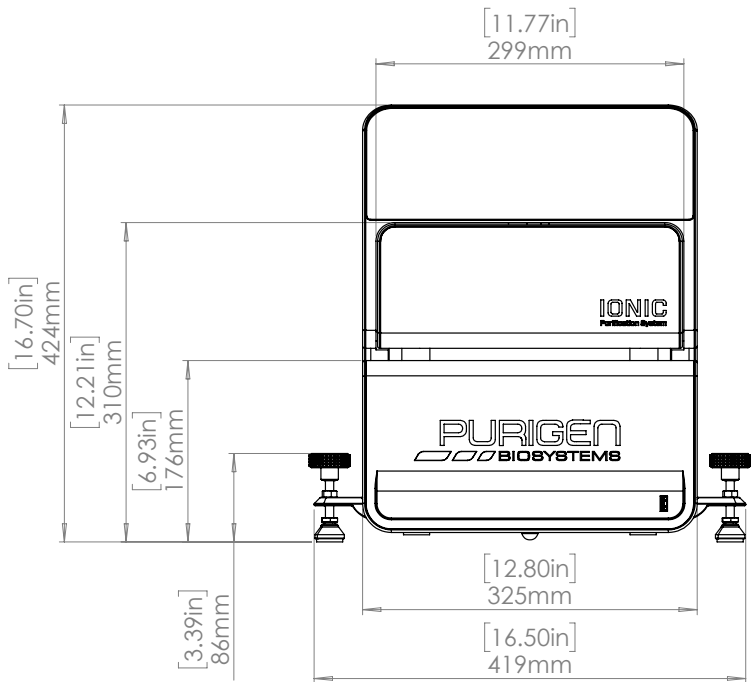


FIGURE 2: Instrument dimensions with leveling system (side)

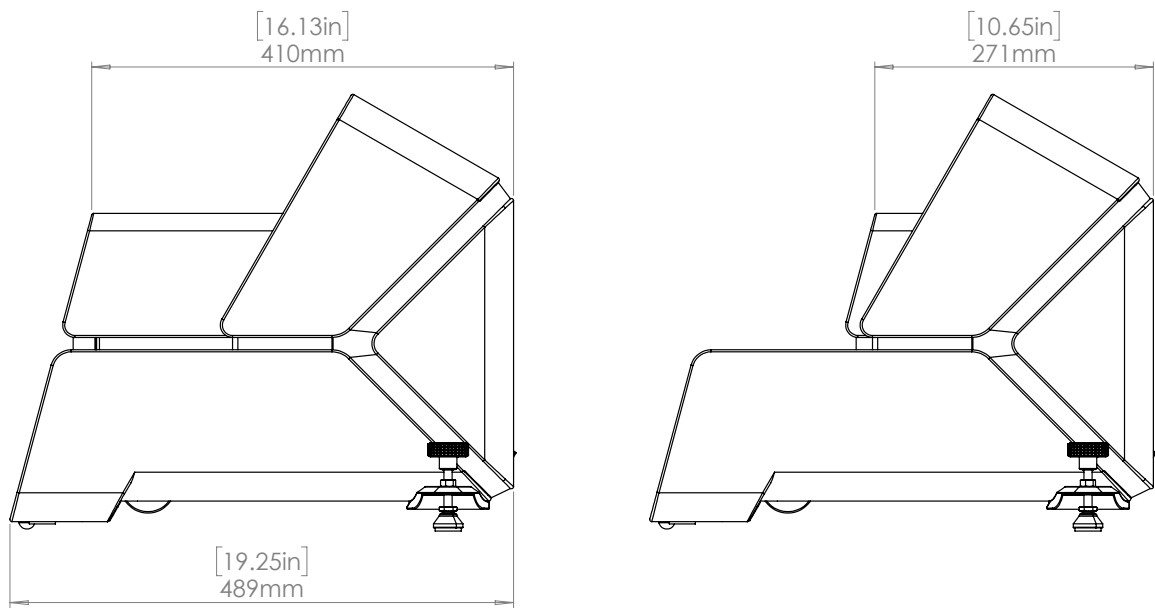

The Ionic system has the following dimensions upon installation:

TABLE 3: Measurement with leveling system

Specification	Value
Height	42.4 cm (16.7 in)
Width	41.9 cm (16.5 in)
Depth	48.9 cm (19.3 in)
Weight	24.5 kg (54 lbs)

Placement Requirements

Select the installation location for your Ionic system prior to scheduled installation. The Ionic system must be positioned to allow access to the power switch and power outlet, for proper ventilation, and for servicing the instrument.

- Make sure that you can reach around the left-side of the instrument to turn on or turn off the power switch on the back panel adjacent to the power cord.
- Position the instrument so that personnel can quickly disconnect the power cord from the outlet.
- The instrument must be accessible from all sides using the following minimum clearance dimensions:

TABLE 4: Measurement with leveling system

Access	Minimum Clearance
Sides	Allow at least 35 cm (14 in.) empty bench space to the left-hand or right-hand side of the instrument for pipettes, pipette tips, reservoir, and keyboard.
Rear	Allow at least 10.16 cm (4.0 in.) behind the instrument.
Top	Allow at least 2 cm (0.79 in.) above the instrument. If the instrument is positioned under a shelf, make sure that the minimum clearance requirement is met.



CAUTION

If you need to relocate the Ionic system, contact your Purigen Biosystems representative. Moving the instrument improperly can impact the alignment and compromise yield and purity.



CAUTION

Heat generated by the instrument is expelled as hot air from the rear of the instrument, which must not be blocked by any object closer than 10 cm to avoid potentially dangerous overheating.

Lab Bench Guidelines


Purigen Biosystems recommends placing the instrument on a fixed lab bench or table capable of supporting an instrument weighing 24.5 kg (54.0 lbs.) and with enough surface area to encompass the instrument and provide adequate workspace for operators. Use of mobile tables or benches with casters should be avoided.



NOTE

Lysate preparation can occur on empty bench space, a separate benchtop, or in a BSL hood.

Vibration Guidelines



CAUTION

The Ionic system is sensitive to vibrations.

Use the following guidelines to minimize vibrations during runs and ensure optimal performance:

- Place the instrument on a sturdy, immobilized surface that is isolated from nearby equipment or designed to suppress the transmission of vibration from nearby equipment.
- Do not place other equipment on the bench that produces vibrations, such as a shaker, vortexer, centrifuge, freezer, or instruments with heavy fans.
- Do not install the instrument near frequently used doors. Opening and closing of doors might induce vibrations.
- Do not install a keyboard tray that hangs below the bench.
- Do not touch the instrument or open the reagent compartment during operation.
- Do not place objects on top of the instrument.

Electrical Requirements

This section lists power specifications and describes electrical requirements for your facility.

Power Specifications

TABLE 5: Instrument parts (front)

Type	Specification
Line voltage	North America / UL120 VAC, 60 Hz, 2.5 A Europe / CE230 V~50 Hz, 1.0 A
Power consumption	300 watts

Recommendations

The following practices are recommended to ensure safe operation of the Purigen Biosystems Ionic system:

- The instrument should be plugged into a receptacle on a dedicated circuit with proper voltage, amperage, and grounding.
- A power line regulator is required if the voltage fluctuates in excess of +/- 10% of the normal voltage.
- A surge protector is recommended in areas with frequent electrical storms or near electrically “noisy” devices such as refrigerators or centrifuges.
- Instrument must be placed less than 2 meters from a power outlet.
- The power outlet must be accessible at all times for safety.

Power Cords

The table below lists the specifications for the power cord included with the instrument.

TABLE 6: Instrument power cord specifications

Specification	North America / UL	Europe / CE
Current rating	10 A	10 A
Voltage rating	125 VAC	250 VAC
Input connection 1	NEMA 5-15P	Various based on country
Input connection 2	IEC 320-C13	IEC 320-C13
Cord length	6.56' (2.00m)	≤ 3m
Cord type	SJT	SVT



WARNING

Never use an extension cord to connect the instrument to a power supply.

Fuses

The Ionic system contains no user-replaceable fuses.

Uninterruptible Power Supply

Purigen Biosystems recommends the use of a user-supplied uninterruptible power supply (UPS). Purigen Biosystems is not responsible for runs affected by interrupted power regardless of whether the instrument is on a UPS. Standard generator-backed power is often not uninterruptible and a brief power outage occurs before power resumes, which interrupts a run.

Purigen Biosystems recommends the use of a UPS rated at least 1500VA/750W with a total capacity of 300Wh to avoid loss of sample and materials during a power outage (APC BX1500M 1500VA Back-UPS Pro, or similar).

Product Certifications and Compliance

The Ionic Purification System is certified to the following standards:

- UL 61010-1, 3rd Edition
- CAN/CSA C22.2 No. 61010-1-12, 3rd Edition
- Subpart B of Part 15 of FCC Rules for Class A digital devices
- ICES-003, Issue 6

European CE Standards

This product fulfills all the relevant provisions of the Low Voltage Directive 2014/35/EC, WEEE Directive and the EMC Directive. The instrument is designed and manufactured in accordance with the following standard(s) noted below.

Title and No. of Directive:	Directive 2014/35/EU on Low Voltage Directive	EN 61010-1:2010
	Directive 2014/30/EU on Electromagnetic Compatibility	EN 61326-1:2013
	WEEE Directive	76/769/EC

Environmental Considerations

The Ionic Purification System is designed to operate within the following environmental conditions:

TABLE 7: Operating environment specifications

Operating Environment	Specification
Temperature	Maintain a lab temperature of 15°C – 25°C. This temperature is the operating temperature of the instrument
Humidity	Maintain a non-condensing relative humidity between 30–70%.

Pre-installation Checklist

- A specific installation site for the system has been selected, is available, and conforms with the requirements specified above
- A primary and secondary operator of the Ionic system has been identified, and contact information has been provided to your Purigen Biosystems Field Service Engineer (FSE)
- Primary Operator Secondary Operator

	Primary Operator	Secondary Operator
NAME		
EMAIL		
PHONE		

The primary and secondary operators will be available on the day of system installation to attend the installation and system training

- Customer-supplied materials listed above are available within arm’s length of the installation site
- Lysate preparation can occur on a benchtop or in a BSL hood

If at the time of scheduled installation the FSE determines that one or more checklist items are incomplete, it may be necessary to re-schedule the installation for another time, and you may be responsible for additional costs.

Your Field Applications Specialist will contact the operator trainees prior to your scheduled installation.

Acknowledgment and Authorization

Please sign and return this form to support@purigenbio.com to confirm that site preparations are complete. Upon receipt, we will return our confirmation of receipt and acknowledgment to arrange installation.

I, _____ agree that I am authorized to approve installation of Ionic Purification System and I agree to have fulfilled all the requirements stated above.

NAME (Signature)	
NAME (Print)	
DATE	
POSITION / TITLE	

For internal Purigen Biosystems support team use only:

I, _____ agree that I have received and approved installation of Ionic system.

NAME (Signature)	
NAME (Print)	
DATE	
POSITION / TITLE	
ADDRESS	Purigen Biosystems, Inc. 5700 Stoneridge Drive, Suite 100 Pleasanton CA 94588 USA